

TRAX

Compression Screw System

Product Information


Breakaway Screw



Headless Screw



Headed Screw

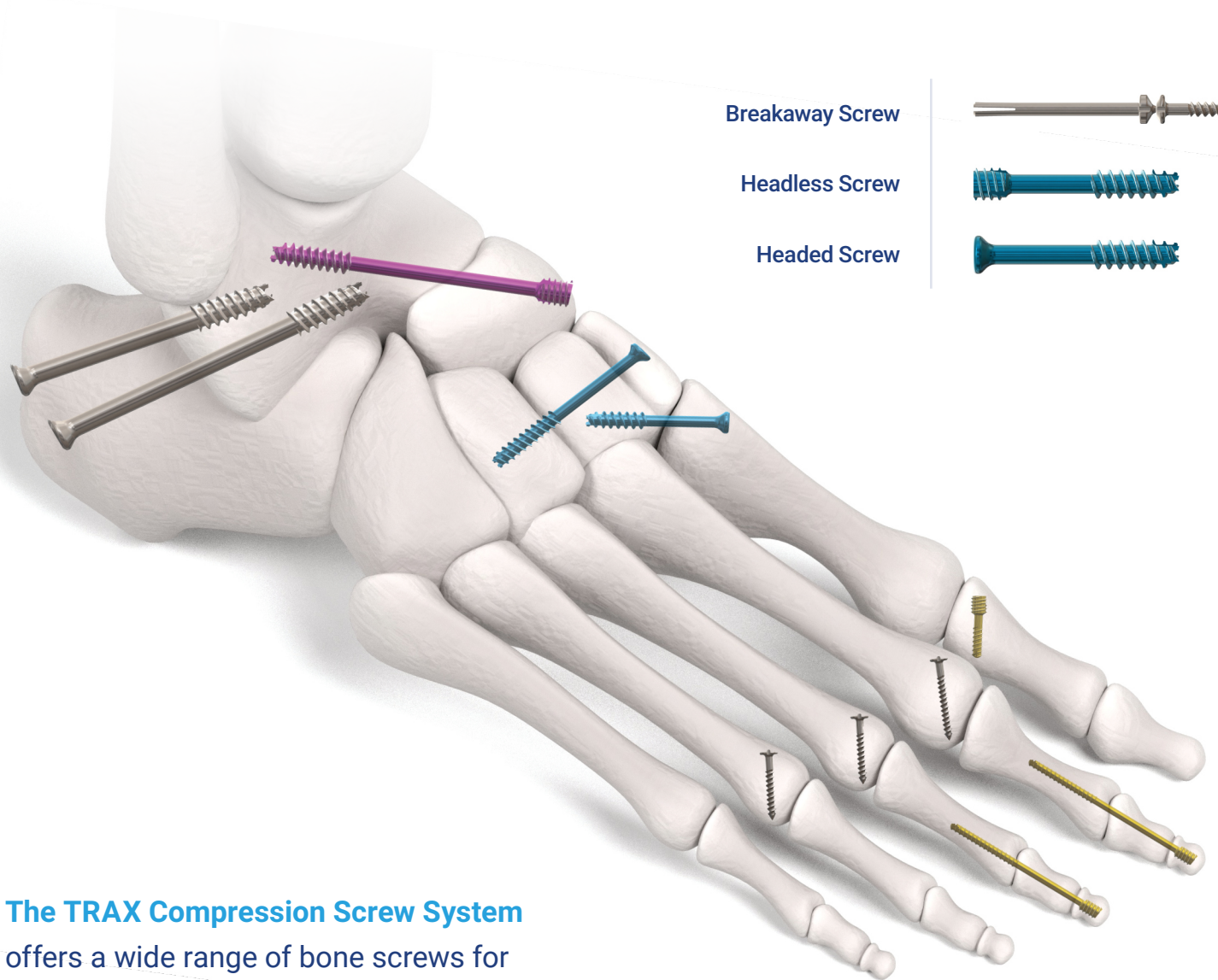

The TRAX Compression Screw System offers a wide range of bone screws for osteotomy, fixation, fusion, and fracture procedures. With several screw and instrument options arranged into three kits, selecting the perfect implant for your patient is straightforward.

TS
TRAX SURGICAL


Features

TRAX Compression Screw

The TRAX Compression Screw System is comprised of headed and non-headed titanium screws. Each canulated, self-tapping screw is available in 2.0mm – 6.5mm diameters, and a choice of lengths that range from 12.0mm – 110.0mm.

Headless Screw Features

Thread design maximizes initial contact and enhances surface area of the screwbone interface.

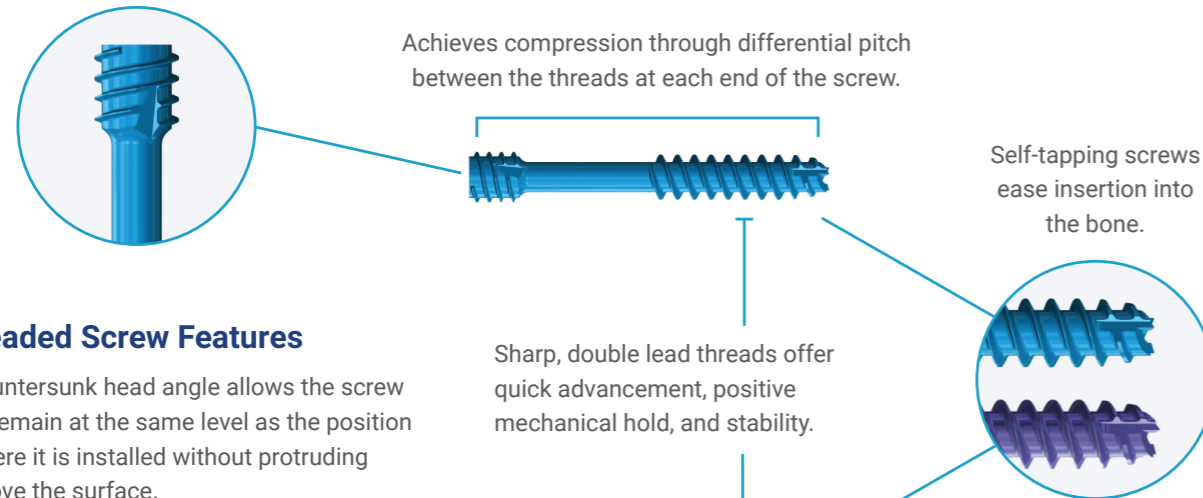


Headed Screw Features

Countersunk head angle allows the screw to remain at the same level as the position where it is installed without protruding above the surface.




The TRAX compression Screw System includes color-coded washers for use with the 3mm, 4mm, 4.5mm, 5.5mm, and 6.5mm headed screws. Washers may be used to optimize compression and minimize the risk of unintentional intrusion of the screw head through cortical bone during screw insertion.

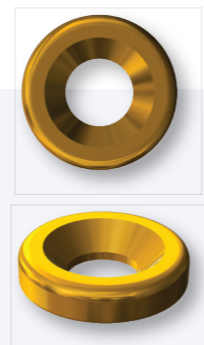

Instrumentation Reference

What do the colors mean? This instrument is compatible with gold and green screws.

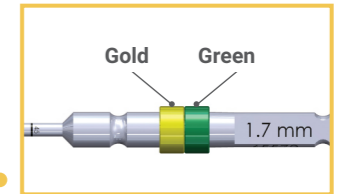

Small Compression Screw Kit REF: ACC-S

GOLD	Headed & Headless 2.0mm Diameter
GREEN	Headed & Headless 2.4mm Diameter
DARK BLUE	Headed & Headless 3.0mm Diameter
LIGHT BLUE	Headed & Headless 4.0mm Diameter

Length (mm)	Thread Length (mm)
* 12	5
12	6
14	7
16	8
18	9
20	10
22	11
24	12
26	12
28	12
30	12
32	16
34	16
36	16
38	16
40	16
42	16
LIGHT BLUE	Headed & Headless 4.0mm Diameter
Length (mm)	Thread Length (mm)
44	16
46	16
48	16
50	16

* Note: applies to headed screws only

Medium Compression Screw Kit REF: ACC-M

ORANGE	Headed & Headless 4.5mm Diameter
MAGENTA	Headed & Headless 5.5mm Diameter

Length (mm)	Thread Length (mm)
* 20	8
* 20	9
20	10
22	11
24	12
26	12
28	12
30	12
32	16
34	16
36	16
38	16
40	16
42	16
44	16
44	32
46	16
46	32
48	16
48	32
50	16
50	32
55	16
55	32
60	16
60	32
65	16
65	32
70	16
70	32
MAGENTA	Headed & Headless 5.5mm Diameter
Length (mm)	Thread Length (mm)
75	16
75	32
80	16
80	32

Large Compression Screw Kit REF: ACC-L

NATURAL	Headed & Headless 6.5mm Diameter
Length (mm)	Thread Length (mm)
40	16
45	16
45	32
50	16
50	32
55	16
55	32
60	16
60	32
65	16
65	32
70	16
70	32
75	16
75	32
80	16
80	32
85	16
85	32
90	16
90	32
95	16
95	32
100	16
100	32
105	16
105	32
110	16
110	32

* Note: applies to magenta, headed screws only

* Note: applies to orange, headed screws only


The TRAX Compression Screw System also includes a set of non-cannulated titanium screws in the small kit. These headed, break-away screws are designed for small bones/fragments, or where small screws are desired for optimum bone reconstruction, repair, osteotomies, and fusions.

**2.0mm Diameter
Break-Away Screws**

Length	Qty/Kit
11mm	4
12mm	4
13mm	4
14mm	4

Threads designed for rapid insertion into cancellous bone, and maintain a firm hold.

Pin breaks away below the screw head to minimize prominence. Lessens risk of soft tissue irritation.

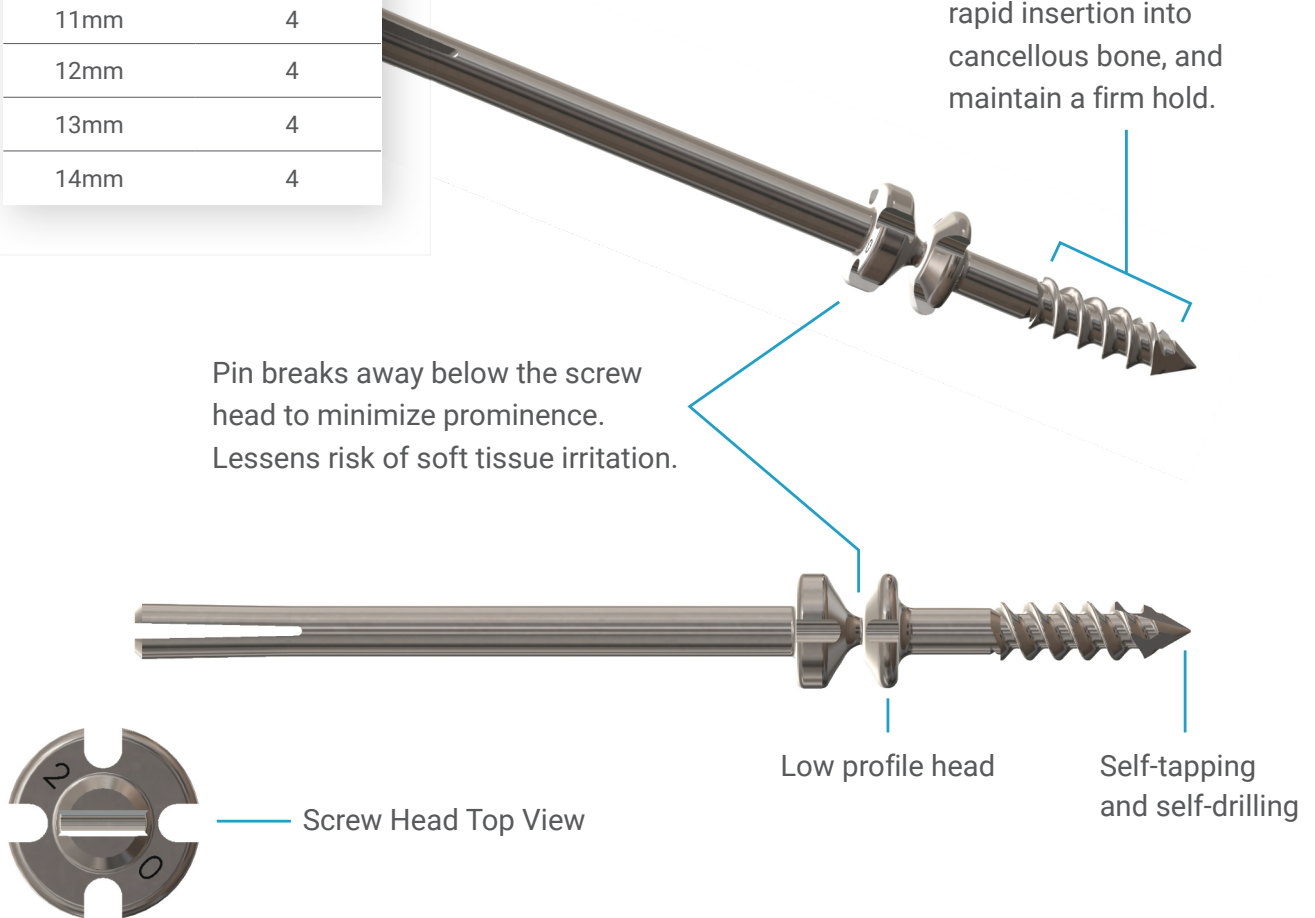

Trax Surgical Inc.

75 Mill Street
Stoughton, MA 02072
Tel: 781-436-4350
www.TraxSurgical.com

P/N 67985 Rev.A 09/2023

